



caring mill by **aura**



Model No. AR60

Revive Ultra Gun

User Manual



**Thank you for purchasing
the Caring Mill by Aura Revive Ultra Gun.**

Prepare to experience ultimate relaxation with adjustable speed levels and interchangeable customized attachments, while advanced technology allows for accelerated muscle recovery, improved blood circulation, and tension relief.

Experience the transformative effects of the Caring Milly by Aura Revive Ultra Gun and embrace a new sense of vitality.

Contents

Warning	4
Package Contents	5
Product Components	6
Interchangeable Attachments	7
Instructions for Use	9
Charging & Battery Information	10
Care & Cleaning Instructions	12
Technical Data/ Specifications	12
Customer Service	13
FCC Compliance Information	13

Save this user manual for future reference



Note

If this User Manual is damaged, or if you no longer have the User Manual in your possession, please contact customer service.

Warning



WARNING

**READ AND FOLLOW THE WARNINGS CAREFULLY
TO REDUCE THE RISK OF EXPLOSION,
FIRE, ELECTRICAL LEAKAGE, OR INJURY.**

- Consult a doctor before use, especially if you have any injuries.
- This product is for adult use only;
keep it away from children and pets.
- Use the device according to the provided instructions.
- Keep the Ultra Gun away from water, heat sources, and ensure the vents remain clear.
- Turn off the appliance when not in use.
- Do not use the gun continuously for more than 30 minutes;
take a five minute break after every 20 minutes of use.

- Avoid using on hard areas such as the head and sensitive areas like the eyes or ears.
- Be cautious when selecting the vibration speed to prevent contusions.


Regularly check the treated area and discontinue use if you


- experience pain or discomfort.
- Use the equipment on dry, clean clothing or skin.
- Keep hair, fingers, and other body parts away from the connection of the interchangeable attachment to prevent pinching.
- Be cautious of prolonged use, as it may cause damage to fascia tissue.

Package Contents

- 1x Ultra Gun
- 1x Charging Dock
- 1x USB Charging Cable
- 10x Interchangeable Attachments
- 1x User Manual

Product Components

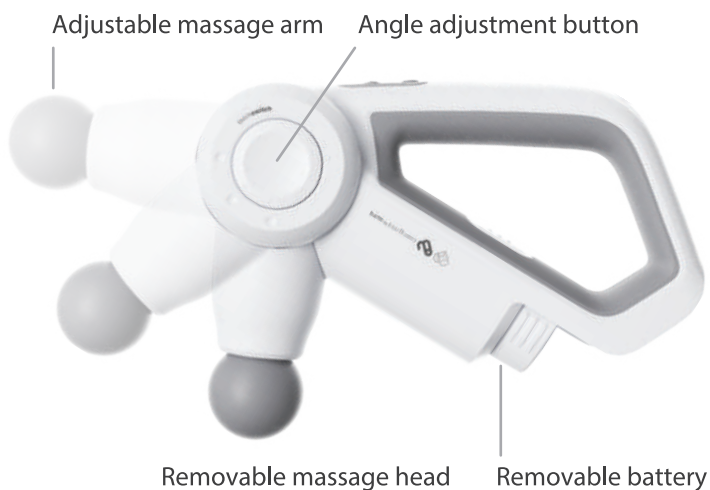

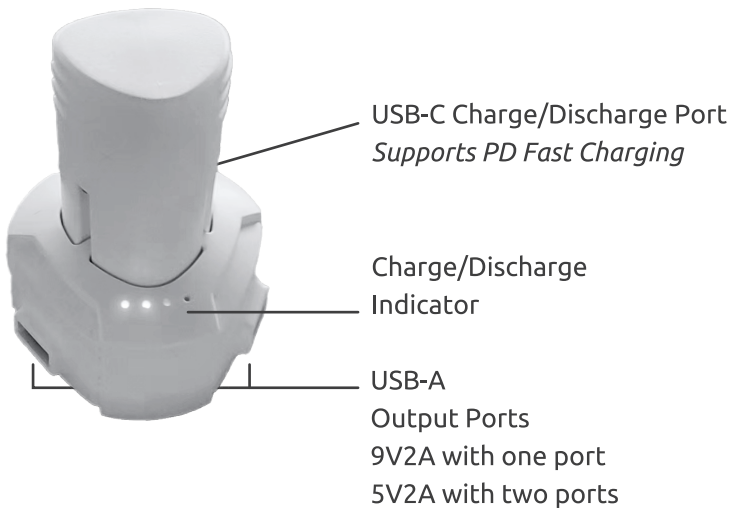

Interchangeable Attachments

Connection Instructions:

- Ensure that the Ultra Gun is fully charged and turned off.
- Choose the desired interchangeable attachment for your session.
- Insert the attachment into the arm of the device.
- Press firmly to securely attach the part, aligning the small projection on the Y-shaped head with the recessed hole in the arm.

Replacement Instructions:

- Power off the Ultra Gun using the power switch.
- Firmly grasp the attachment currently on the body of the device.
- Pull the attachment out firmly to remove it from the arm.
- Insert the new attachment that you wish to use.
- Press firmly to ensure the part is securely attached and not easily removable.



Whishbone

The wide wishbone design provides a strong massage around joints and hard angles, reducing hand fatigue and awkward angles.



Kayak

The wide curved head glides gently over joints and hard angles which reduces hand fatigue and awkward angles.



Precision point

The small point on this massage head allows you to precisely pinpoint problem areas.



Flat

For large muscle groups and areas that need a little more help. Simulates a deep-tissue massage.



Round

Improve recovery and reduce muscle soreness after workouts with gentle oblique impact.



Big Soft

This attachment is designed specifically for low-impact treatment around tender or bony areas.



Textured Flat

The textured surface helps loosen tight muscles. For large muscle areas that need extra attention.



Shovel

For back muscles and areas surrounding the spinal column. Great for flushing out toxins.



Spinal

Loosen up tough back tissue and improve flexibility with the spinal massage head.



Bullet

Designed specifically to target problem areas and release muscle knots.

Instructions for Use

- 1 To power on the Ultra Gun, press and hold the power switch. Once the power is on, press the speed control button to initiate the device's operation.
- 2 Adjust the speed of the massage gun by pressing the speed control button in the following sequence:
1 (LED green) - 2 (LED blue) - 3 (LED purple) - 4 (LED cyan) - 5 (LED red).



- 3 If you need to pause the session while operating between gears 1 and 5, press the power switch. The next time you start the Ultra Gun, it will automatically resume at the previously set speed.
- 4 When you have completed use, press the power switch to turn off the device and stop its operation.

Charging & Battery Information

When the Ultra Gun is powered on, the power indicator light provides a visual representation of the battery level as follows:


- Green: Indicates battery power between 70% and 100%.
- Blue: Indicates battery power between 40% and 70%.
- Red: Indicates battery power between 10% and 40%

When the power indicator light flashes red, it signifies that the battery is below 10% and requires recharging.



Instructions for Charging the Battery Pack

- 1 Locate the battery eject button on the Ultra Gun.
- 2 Press the battery eject button to automatically eject the battery from its slot.

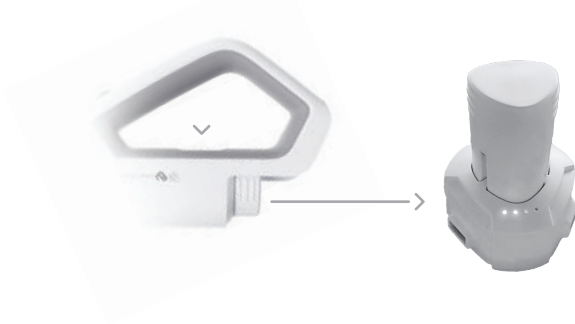

- 3 Remove the battery from the massage gun.
- 4 Connect the data charging cable to the charging cradle's Type C port.
- 5 Insert the battery pack into the charging cradle.
- 6 Ensure a secure connection between the charging cradle and the battery pack.
- 7 Allow the battery pack to charge fully, which typically takes around 3 hours.
- 8 Once the battery is fully charged, align it correctly and insert it directly into the battery slot.
- 9 Ensure a secure fit by listening for the "click" sound, indicating that the battery is properly installed.

Note: For optimal device performance, ensure that the battery pack is fully charged prior to use.

Charging Dock LED Indicator

Power Level	LED 1	LED 2	LED 3	LED 4
Power <30%	FLASH	FLASH	FLASH	FLASH
Power >30%	ON	FLASH	FLASH	FLASH
Power >60%	ON	ON	FLASH	FLASH
Power >80%	ON	ON	ON	FLASH
Power 100%	ON	ON	ON	ON

Care & Cleaning Instructions

Care, Cleaning, and Storage

- Use a soft, clean towel or cloth to wipe away any dust or stains.
- Store or transport the gun in the provided carrying case.

Technical Data / Specifications

Name: Caring Mill by Aura Revive Ultra Gun

Model No AR60

Motor: 12V Brushless Motor

Noise: 60db

Stroke: 16mm

Motor Power: 250W

Weight: 1.59 lbs.

Product Size: 9.05in*5.19in

Rated Speed: 1200-2400 rpm

Battery Capacity: 2000mAh

Charging Method: CC-CV

Rated Input: 5V/2A or 9V/2A

Output Voltage for Battery: $12.6 \pm 0.3V$

Output Voltage for Dock: 5V/2A or 9V/2A

Operating Temperature: 32~131 F

Manufactured for: Aura Medical LLC

2 Skillman St.

Suite 512

Brooklyn, NY 11205

Caring Mill by Aura

Customer Service

For support or inquiries please contact us at:
support@aura-medical.com

FCC Compliance Information



The Caring Mill by Aura Revive Ultra Gun AR60
FCC Compliance Statement
This device complies with part 15 of the FCC Rules.

